

THE + FEATURES OF OSTEO SAFE®

► Calibrated monitoring of surgical procedure. Enhanced comfort and safety for the patient.

- Atraumatic impaction.
- Mini-invasive bone preparation.
- Single hand-held instrument, meaning improved visibility during surgery.
- Monitoring and control of the impaction (speed, force, axis).
- Reproducibility and accuracy of the procedure.
- Simplified protocol using 4 osteotomes for placement of Axiom® REG or PX implants (intraoperative time-saving).



*For more information on Osteo Safe®,
please contact your Anthogyr local representative.*

Manufacturer: Anthogyr. Please carefully read the instructions given in the notices and user guides.

ANTHOGYR SAS
2 237, Avenue André Lasquin
74700 Sallanches - France
Tél. +33 (0)4 50 58 02 37
Fax +33 (0)4 50 93 78 60
www.anthogyr.com

Anthogyr
PRIME MOVER IN IMPLANTOLOGY

SAFETY - REPRODUCIBILITY - ACCURACY

OSTEO SAFE®

NEW!



Performing calibrated
osteotomy safely

Anthogyr
PRIME MOVER IN IMPLANTOLOGY

C165_GB-2015-02

Photo Credits: Anthogyr – All rights reserved – Non contractual photos.

- ▶ Anthogyr has developed an innovative solution for easy and calibrated crestal impaction osteotomy: **Osteo Safe®**.

UNIQUE!

Precalibrated instrument for automatic impaction. Designed to be connected to a standardised universal micromotor.



The Osteo Safe® can be used **for all indications of implant site preparation and bone remodeling** in the context of a vertical and lateral bone augmentation procedure :

- Crestal sinus lift.
- Bone compaction in order to obtain sufficient primary stability for placement of Axiom® REG/PX implants (low bone density in the upper-jaw: *bone type: III or IV).



In order to provide an appropriate response to daily practice, the Osteo Safe® and its instruments are supplied **in a sterilisable surgical cassette** containing:

- 1 Osteo Safe® impactor,
- 4 straight osteotomes (insert tips),
- 4 bayonet osteotomes (insert tips),
- 4 **concave-shaped** insert tips (A),
- 4 **convex-shaped** insert tips (B),
- 1 manual screwdriver and its accessories (mandrel and mandrel extension).
- **The conception of the Osteo Safe® system enables to customize the kit with appropriate modules.**

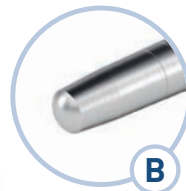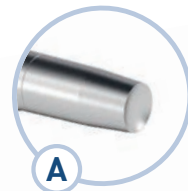

An update for Implanteo® will soon be available at www.anthogyr.com



*according to LEKHOLM & ZARB Classification (1985).